

|                                                                                                                            |                                                                                         |                                                      |                                                   |
|----------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------|------------------------------------------------------|---------------------------------------------------|
|  <b>DIPLOMAT<br/>DENTAL<br/>SOLUTIONS</b> | <b>Diplomat Dental s.r.o ., Vrbovská cesta 17, 921 01 Piešťany, Slovenská republika</b> |                                                      |                                                   |
|                                                                                                                            | <b>Change description :</b>                                                             | <b>Change is valid for :</b><br>Model PRO, Model ONE | <b>Technical information</b><br><b>No.: 11_26</b> |
|                                                                                                                            | Modification of the handle mounting design on SIRIUS operating lights                   | <b>Change is valid from :</b><br>08/10/2026          | <b>Revision:</b>                                  |


## Modification of the handle mounting design on SIRIUS operating lights

Dear Customer, we would like to inform you about a structural change regarding the **SIRIUS** family of operating lights (**ESSENTIAL, PLUS, POWER, and ULTRA models**), intended for the **MODEL PRO** and **MODEL ONE** dental units.

This change focuses on modifying the mounting mechanism for both the **FIX HANDLES (ME348000.071)** and the **REMOVABLE HANDLES (ME348000.072)**.

### Details of the technical change:

To simplify assembly and ensure the correct seating of the handles, a new structural design has been developed. New positioning projections have been added to the aluminum castings of the light body. At the same time, corresponding guiding grooves have been created on the plastic moldings that are part of the handles. This eliminates the previous design, which was more demanding for precise assembly.

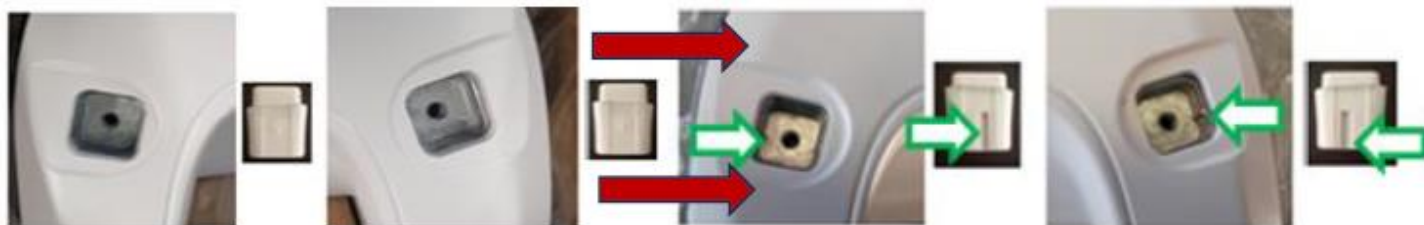

### Questions and answers:

**Question:** Does this structural modification change the way removable handles are cleaned or sterilized?

**Answer:** No, the change is purely structural and only affects the internal mounting mechanism. The maintenance, cleaning, and sterilization processes for removable handles remain completely unchanged.

**Question:** I have an older model of the SIRIUS light in my clinic without these positioning projections. Do I need to replace the handles?

**Answer:** No, if your current light and handles are functioning correctly, no replacement is necessary. This modification applies only to new deliveries from the factory to simplify installation and increase the overall precision of the part fitting.

**Question:** Are the new handles compatible with older versions of the SIRIUS lights and vice versa?

**Answer:** No, due to the structural changes (addition of positioning projections and guiding grooves), the new handles are not compatible with older lights. Likewise, older handles cannot be mounted on lights with the

| Diplomat Dental s.r.o ., Vrbovská cesta 17, 921 01 Piešťany, Slovenská republika                                           |                                                                       |                                                      |                                             |
|----------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------|------------------------------------------------------|---------------------------------------------|
|  <b>DIPLOMAT<br/>DENTAL<br/>SOLUTIONS</b> | <b>Change description :</b>                                           | <b>Change is valid for :</b><br>Model PRO, Model ONE | <b>Technical information<br/>No.: 11_26</b> |
|                                                                                                                            | Modification of the handle mounting design on SIRIUS operating lights | <b>Change is valid from :</b><br>08/10/2026          | <b>Revision:</b>                            |

new structural design. When ordering spare parts, it is therefore important to specify which version of the light you have on your dental unit.

If you have any further questions regarding this technical information or any other inquiries, please contact us via our official ticket portal at [ticket.diplomat-service.com](https://ticket.diplomat-service.com).

Thank you,

Your Diplomat Dental Solutions team.